LiFE

CARDIO

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PLEASE READ AND SAVE THESE IMPORTANT SAFETY INSTRUCTIONS.

Read the operating instructions carefully before putting the appliance into operation and keep the instructions, the receipt and, if possible, the box with the internal packing. If you give this device to other people, please also pass on the operating instructions.

USER MANUAL SYMBOLS

Important information for your safety is specially marked. It is essential to comply with these instructions in order to avoid accidents and prevent damage to the machine:

 **WARNING:** This sign warns you of dangers to your health and indicates possible injury risks.

 **CAUTION:** This sign refers to possible hazards to the machine or other objects.

 **NOTE:** This sign highlights tips and information

SAFETY INSTRUCTIONS

1. The appliance is designed exclusively for private use and for the envisaged purpose.
2. This appliance is not for professional use.
3. Do not use the scale outdoors.
4. Keep the appliance away from sources of heat, direct sunlight, humidity (never dip it into any liquid) and sharp edges.
5. Do not try to repair the appliance on your own. Always contact a qualified technician.
6. Do not overload the scale by exceeding the maximum weight of 180kg as the sensors of the appliance may be damaged.
7. In order to ensure your children's safety, please keep all packaging (plastic bags, boxes, polystyrene etc.) out of their reach.
8.  **WARNING:** Do not allow small children to play with the packaging as there is a danger of suffocation!
9. Children should be supervised in order to ensure that they do not play with the appliance.
10. When the display shows "Lo" it is indicated that the power of the batteries is very low. Please replace the batteries. Low battery power or defective


- batteries may cause inaccurate results or other errors.
11. Use the appliance only on a hard, flat and dry surface.
 12. Do not place the scale under wet, hot & cold conditions.
 13. Clean the appliance with a damp cloth and a mild detergent.
 14. Avoid excessive impact or vibration on the appliance.
 15. If you do not intend to use the appliance for a long time, remove the batteries before storage to avoid a battery fluid leakage.

SPECIFICATIONS

1. Glass digital bathroom scale for analysis of weight, body fat content, water content, muscle mass content, bone weigh content, calorie requirements per day (BMR) and BMI value (Body Mass Index)
2. Large LCD display (76x40mm) with large and legible white LED digits
3. Memory for up to 10 different users
4. Measurement range: 5-180kg
5. Maximum loading: 180kg
6. Accuracy: 100g
7. Step-on-automatic
8. Large high-safety glass surface 30x30cm, made of 6mm thick tempered glass
9. Automatic switch off, prolongs battery life
10. Overload indicator
11. Low battery indicator
12. Rubber feet for maximum stability
13. Battery powered

PARTS IDENTIFICATION

1. Display
2. Weighing platform
3. SET button
4. UP button ▲
5. DOWN button ▼

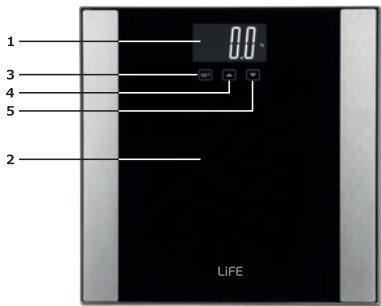LCD display

- 1. Selection of the type of personal physical condition of the user (athletic / normal)
- 2. Indication of muscle mass (10% ~ 80%)
- 3. Man / Woman
- 4. Bone mass content (0.5Kg ~ 8.0Kg)
- 5. Age (range 10 - 80 years)
- 6. Indication of body fat (5% ~ 50%)
- 7. Indication of daily need in calories
- 8. Indication of body hydration (10% ~ 80%)
- 9. User memory location
- 10. Body mass index (BMI)
- 11. Indication of Underfat / Healthy / Overfat / Obese
- 12. Height
- 13. Units of measurement

INSERTING THE BATTERIES

- 1. Open the battery compartment lid on the bottom of the scale.
- 2. Insert the batteries into the battery compartment with the correct polarity.


⚠ WARNING:

- Do not expose the batteries to high temperatures or direct sunlight.
- Keep batteries away from children. Batteries may be swallowed, which can be life-threatening. If a battery has been swallowed, seek medical help immediately.
- If the device is not used for a prolonged period of time, please remove the batteries in order to prevent leakage of battery acid.

USING THE APPLIANCE

SELECTION OF WEIGHT MEASUREMENT UNIT

Select the desired measurement unit in kg, lb. or st, by pressing the “UNIT” button on the bottom side of the scale. The default setting is kg.

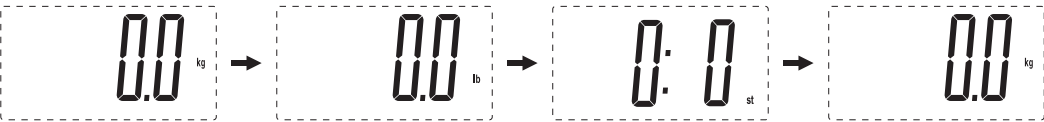
PERSONAL INFORMATION

You must enter the following user parameters to measure body fat percentage and other data.


- 1. Place the appliance on a hard, flat, dry surface. Avoid placing it on the carpet or on soft surfaces.
- 2. Turn on the scale and press the “SET” button.

3. Select the user memory spare using the  and  confirm the setting by pressing the "SET" button.

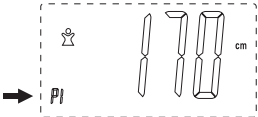4. Select the gender ("Male" or "Female") and your fitness level (athletic or normal) using the  and  buttons and confirm the setting by pressing the "SET" button.

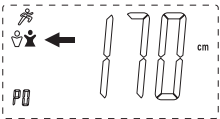i NOTE:

- A person who is athletic can be defined as a person:
 - Aged 18 or above
 - Takes part in more than 6 hours of intense aerobic exercise weekly
 - Has a resting heartbeat below 60 beats per minute.
- If you select a sport activity, all measurements are affected.
- We do not recommend using the sports activity feature unless you meet the above requirements.

5. Select your height using the  and  confirm the setting by pressing the "SET" button.

6. Select your age using the  and  confirm the setting by pressing the "SET" button.

WEIGHING


1. Place the appliance on a hard, flat, dry surface. Avoid placing it on the carpet or on soft surfaces.
2. Switch on the scales by pressing lightly with your foot or toes on the surface and wait a few moments. "0.0" will appear on the display.

3. Stand on the scale with both feet, the indicators will flash for a while and then your weight will appear on the screen.


4. Then, the four “0” disappear one by one from left to right.

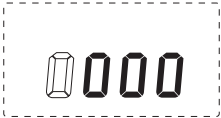
5. After a few seconds, the following results will appear:

- your weight percentage
- body fat percentage in % (FAT symbol)
- body hydration percentage in % (TBW symbol)
- muscle mass percentage in % (♂ symbol)
- bone mass content percentage in% (♀ symbol)
- calorie requirements per day (KCAL)
- body mass index (BMI)

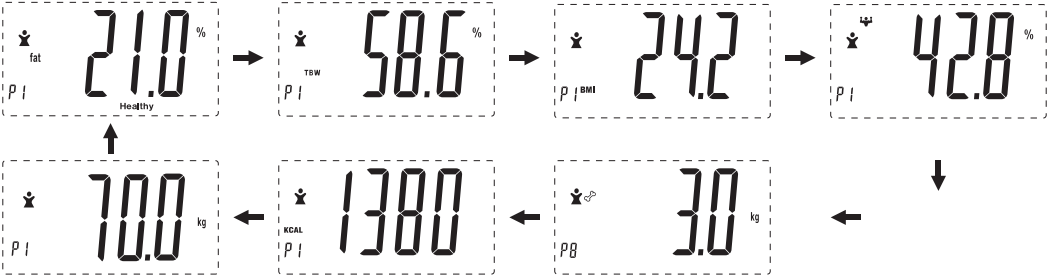
To be able to see all the results, the scale will display them twice before turning off. At the end, it will display the user information (height and age) for confirmation. Then the scale will turn off automatically.

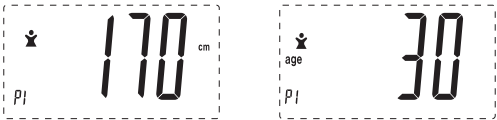AUTOMATIC USER IDENTIFICATION

If the change of your weight compared to the last time measurement is within ± 2 kg, the scale recognizes the user and automatically weighs and analyzes the measurements.

i NOTE: If 2 or more users have the same weight, use the  and  keys to select the correct user.

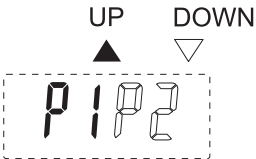i NOTE: If you have not saved a user or there is no automatic user identification, then the screen will display your weight and then turn off automatically after a few seconds.

EVALUATION OF RESULTS – PARAMETERS REFERENCE

BODY FAT (indicative values)

Body fat analysis (Unit %)								
Age	Female				Male			
	Underfat	Healthy	Overfat	Obese	Underfat	Healthy	Overfat	Obese
10-39	<21	21-33	33-39	>39	<8	8-20	20-25	>25
40-59	<23	23-34	34-40	>40	<11	11-22	22-28	>28
60-80	<24	24-36	36-42	>42	<13	13-25	25-30	>30

BODY HYDRATION (indicative values)

Water is a key component of the human body, but the percentage of body water varies depending on age and gender. However, it is important to maintain a body water content above 50%.

Body water analysis (Unit %)						
Age	Female			Male		
	Underfat	Healthy	Overfat	Underfat	Healthy	Overfat
10-15	<57	57-67	>67	<58	58-72	>72
16-30	<47	47-57	>57	<53	53-67	>67
31-60	<42	42-52	>52	<47	47-61	>61
61-80	<37	37-47	>47	<42	42-56	>56

BODY MASS INDEX (BMI) (indicative values)

The Body Mass Index is widely used to assess whether a person has a healthy weight for his height, as, BMI is an index based on the weight and height of individuals. With the BMI index, a person can assess whether he is at risk for health problems, such as heart disease, diabetes, etc. The table below gives an indicative report on the assessment of health status using BMI.

BMI parameter (reference to WHO)				
	Under weight	Normal	Over weight	Obese
BMI	<18.5	18.5~25	25~30	>30

BODY MUSCLE MASS PERCENTAGE (indicative values)

Body muscle analysis (Unit %)		
Gender	Female	Male
Normal	>34%	>40%

BODY BONE MASS PERCENTAGE (indicative values)

Body bone analysis (Unit Kg)						
Gender	Female			Male		
Weight	<45Kg	45-60Kg	>60Kg	<60Kg	60-75Kg	>75Kg
Bone mass	1.8Kg	2.2Kg	2.5Kg	2.5Kg	2.9Kg	3.2Kg

DAILY NEED CALORIES (indicative prices)

Calories analysis (Unit:KJ)					
	Age	<17	18~29	30~50	>50
Female	<50Kg	1265	1180	1085	1035
	50~65Kg	1392	1298	1194	1139
	65~80Kg	1771	1652	1519	1449
	80~95Kg	2151	2006	1845	1760
	>95Kg	2530	2360	2170	2070
Male	<50Kg	1350	1200	1115	1075
	50~65Kg	1485	1320	1227	1183
	65~80Kg	1890	1680	1561	1505
	80~95Kg	2295	2040	1896	1828
	>95Kg	2700	2400	2230	2150

CLEANING

- 1. ⚠ **CAUTION:** Do not immerse the scale in water: the electronic components can be damaged.
- 2. Do not use a wire brush or any abrasive items.
- 3. Do not use any acidic or abrasive detergents.
- 4. Clean the appliance with a slightly moist cloth. Make sure that no humidity or waters seeps into the device!
- 5. Use a commercially sold detergent or soap for cleaning.
- 6. Dry the device well after cleaning with a dry and soft cloth.

STORAGE

- 1. Always store the appliance away from children.
- 2. Always store the appliance in a safe and clean place. For the best protection of the device, store it in its original packaging.

DISPOSAL OF APPLIANCE AND ACCESSORIES

 This symbol appearing on the product, its accessories or accompanying manuals indicates that the product and its electronic accessories should not be disposed of with other household trash.

In order to avoid possible harmful effects on the environment or health due to uncontrolled waste disposal, please separate these products from other types of waste and recycle them. Household users should contact either the store where they purchased this product or their local services for details of where and how they can return these products for environmentally safe recycling.

DISPOSAL OF BATTERIES

This marking on the battery, manual or packaging indicates that the batteries of this product should not be disposed of with other household waste at the end of their life cycle. ion of the aquifer and subsequent contamination of the food chain causing health problems.